


OTICON | **Real**
OTICON | **More**
OTICON | **Own**
OTICON | **Play PX**
OTICON | **Zircon**

Easy, hands-free communication

Delivered via Bluetooth® Low Energy technology, your patient can now participate in hands-free calls on select iPhone® or iPad®* using their Oticon hearing aids as the headset.

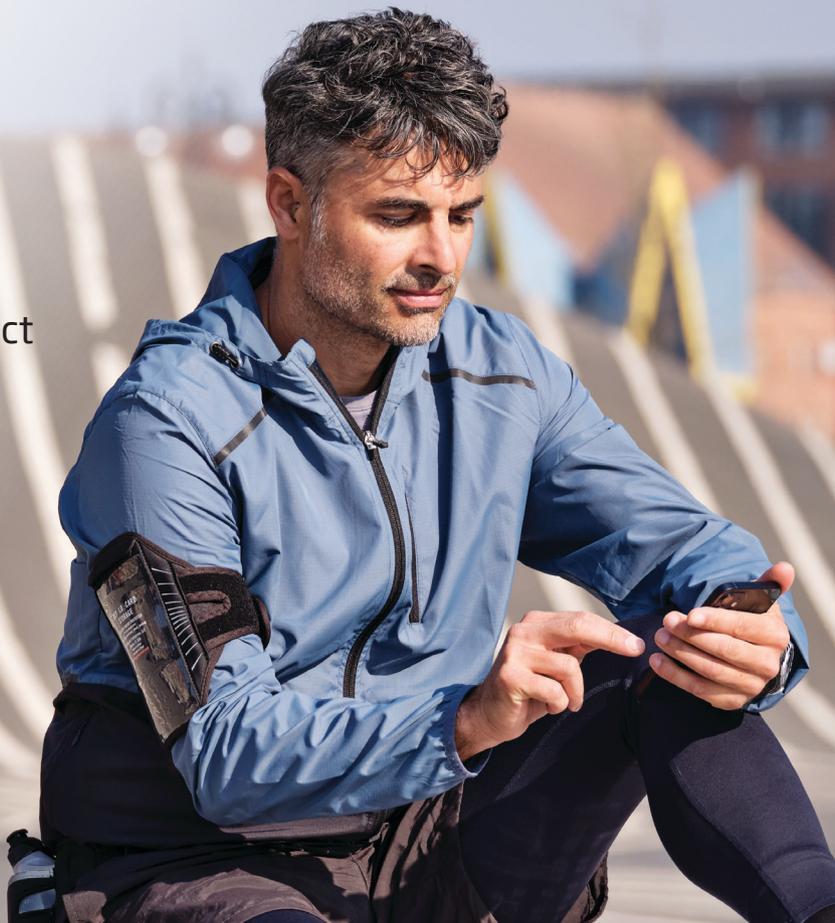

Answer via
push button**
for even easier
call handling



More convenient than ever

Going hands-free lets your patients enjoy the simple things – like calling their loved ones on-the-go, or taking important business calls in the car.

Most importantly, high-quality sound and microphone pick-up allow them to stay focused, multi-task, and let their hearing aids take care of the rest.

- Participate in phone calls and video calls with high-quality, two-way audio
- More convenience when taking remote calls
- Enables better multi-tasking for your patient

Made for
iPhone | iPad

* Hands-free communication is only available in styles with 2.4 GHz connectivity. Hands-free communication is available with iPhone 11 or later running iOS 15.2 or later, and iPad running iPadOS 15.2 or later

** For hearing aids with single or double push button and 2.4 GHz

Apple, the Apple logo, iPhone, iPad, and iPod touch are trademarks of Apple Inc., registered in the U.S. and other countries. App Store is a service mark of Apple Inc.

oticon
life-changing technology

How it works

1. The wearer makes or receives a call on their iPhone or iPad
2. In-call audio is sent directly to their hearing aids
3. The hearing aid microphones capture the wearer's voice
4. The wearer's voice is sent directly back to their device



Built on the latest Bluetooth® Low Energy technology

Oticon's hands-free communication uses less power than classic Bluetooth – enabling your patients to enjoy longer battery life when making and receiving calls.




Powered by the Polaris R™ and Polaris™ platforms

Hands-free communication is available for Oticon Real™, Oticon Own™*, Oticon More™, Oticon Zircon, and Oticon Play PX. Oticon More requires 1.4.0 firmware update for push button control; Oticon Play PX and Oticon Zircon require 1.2.0 firmware update for push button control.

See which hearing aids and devices are compatible here
[oticon.com/support/compatibility](https://www.oticon.com/support/compatibility)

* Only available in styles with 2.4 GHz technology.

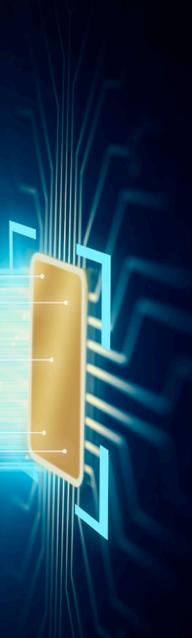

22-356700_15500-11114/03.23