

Jet PX 1

Speech Understanding

OpenSound Navigator™	•
Balancing power effect	40%
Max. noise removal complex/simple	6 dB / 0 dB
Speech Guard™	•
Frequency lowering	Speech Rescue™

Sound quality

Fitting Bandwidth ¹	8 kHz
Power Bass (streaming)	•
Processing Channels	48

Listening comfort

Feedback Management	SuperShield & Feedback shield
Transient Noise Management	On/Off
Wind Noise Management	•

Personalization & optimized fitting

Fitting Bands	14
Multiple Directionality options	•
Adaptation Management	•
Fitting Formulas	NAL-NL1/ NAL-NL2, DSL v5

Connecting to the world

Hands-free communication ²	•
Direct streaming ³	•
Oticon Companion app	•
ConnectClip	•
EduMic	•
Remote Control 3.0	•
TV Adapter 3.0	•
Phone Adapter 2.0	•
Tinnitus SoundSupport™	•
CROS/BiCROS support	•



• Indicates feature included

¹Bandwidth accessible for gain adjustments during fitting


²Hands-free communication is available on select devices

³From iPhone®, iPad®, Mac® and select Android™ devices

Fitting ranges - miniRITE

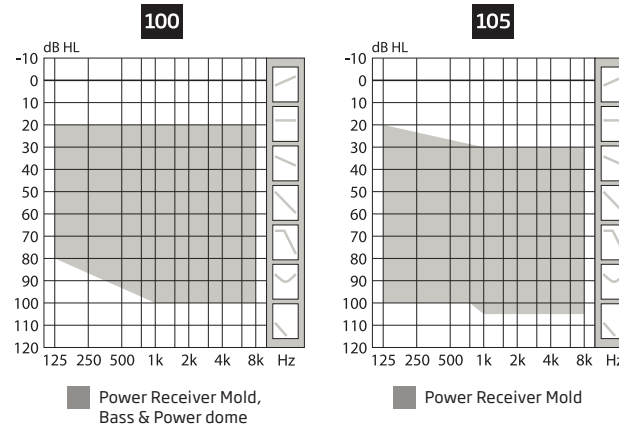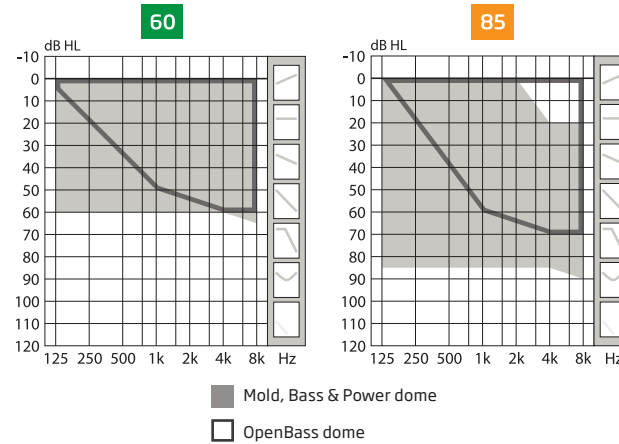

Fitting ranges - miniBTE

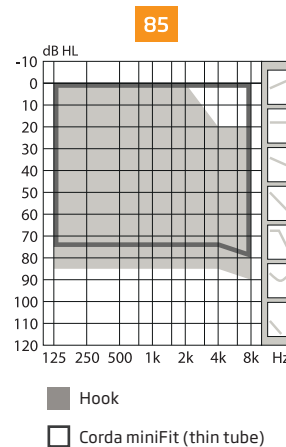

OTICON | Jet PX

Product comparison



oticon
life-changing technology

Oticon Jet PX

The world's first essential level hearing aid with Deep Neural Network (DNN) technology

Oticon Jet PX is powered by the Polaris™ platform and redefines the essential hearing care category.

- Features second generation OpenSound Navigator™ with AI-enabled DNN technology.
- 2.4 GHz Bluetooth® Low Energy technology for wireless connectivity
- Hands free Communication¹
- Direct Streaming²

Oticon Jet PX miniRITE R - is powered by a rechargeable lithium-ion battery, and offers a discreet design with LED-light to make handling easy. The style features a telecoil, and a double push-button.

Oticon Jet PX miniRITE T - is powered by a disposable 312 zinc-air battery, and offers a discreet design with LED-light to make handling easy. The style features a telecoil, and a double push-button.

Oticon Jet PX miniBTE R - is powered by a rechargeable lithium-ion battery, and offers a discreet design with LED-light to make handling easy. The style features telecoil, and a single push-button.

Oticon Jet PX miniBTE T - is powered by a disposable 312 zinc-air battery, and offers a discreet design with LED-light to make handling easy. The style features telecoil and a single push-button.



Apple, the Apple logo, iPhone, iPad, Mac and the Mac logo are trademarks of Apple Inc., registered in the U.S. and other countries. Use of the Made for Apple badge means that an accessory has been designed to connect specifically to the Apple product(s) identified in the badge, and has been certified by the developer to meet Apple performance standards. Apple is not responsible for the operation of this device or its compliance with safety and regulatory standards. Android™ is a trademark of Google LLC. The Bluetooth® word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. and any use of such marks by Demant is under license. Other trademarks and trade names are those of their respective owners.

¹Hands-free communication is available on select devices

²From iPhone, iPad, Mac and select Android devices


For information on compatibility, please visit oticon.com/support/compatibility

oticon.com



miniFit receiver options:

Select between receivers with four different output performances.

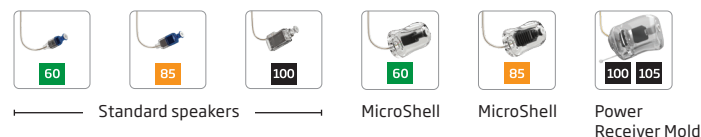


miniFit earpieces:

miniFit earpieces include standard domes and custom molds. LiteTip is available in acrylic or OtoTherm™ material. MicroMold is available in acrylic, soft silicone, or OtoTherm material. MicroShell and Power Receiver Mold are only available in acrylic.



Hook and thin tube for miniBTE:

Select between miniBTE hook with one of our custom molds or Corda miniFit thin tubes with one of our miniFit earpieces.

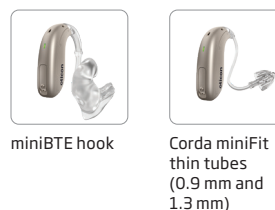

oticon
life-changing technology