

## PRODUCT INFORMATION OTICON CONNECTLINE PHONE 2.0

*The **Oticon ConnectLine Phone Adapter 2.0** is an interface between a conventional landline/home phone installation and the Oticon Streamer, enabling wireless, hands-free phone conversations. The user's Oticon hearing instruments effectively become a wireless headset.*

*As part of the ConnectLine system, Phone Adapter 2.0 works seamlessly and simultaneously with other ConnectLine solutions such as the TV Adapter and the ConnectLine Microphone.*

### **A fully-integrated system**

With ConnectLine, the user is in control of the system. Advanced features such as AutoPriority and AutoConnect ensure that the ConnectLine system works intuitively. For example, users can hear and answer a phone call while watching TV.

### **Landline home phone**

Phone Adapter 2.0 is installed in parallel to an existing landline home phone installation.

Once a user is in range of the Phone Adapter, the Streamer is automatically connected and ready to receive calls (AutoConnect). Incoming phone calls are easily answered with a single button press on the Streamer at a distance of up to 30 meters from ConnectLine Phone 2.0.

### **AutoPriority**

If a user is listening to the TV using the ConnectLine TV solution (or other audio sources), the AutoPriority feature ensures that an incoming call can be heard and answered.

### **AutoResume**

When a call is ended, Streamer automatically returns to its previous state. For example, back to TV input if a call is taken while watching television.

For more information please refer to:  
[www.oticonusa.com/connectline](http://www.oticonusa.com/connectline)

### **Office phone**

Phone Adapter 2.0 can be used in an office environment to create a wireless headset solution for digital/PABX-based phones. The additional components needed for these applications are available in a ConnectLine office phone kit.

For more information please refer to:  
[www.oticonusa.com/connectline](http://www.oticonusa.com/connectline)




## PRODUCT OVERVIEW

# Oticon | ConnectLine

### Included

- RJ11 - RJ11 cable for connection to landline home phone
- Power supply
- Instructions For Use

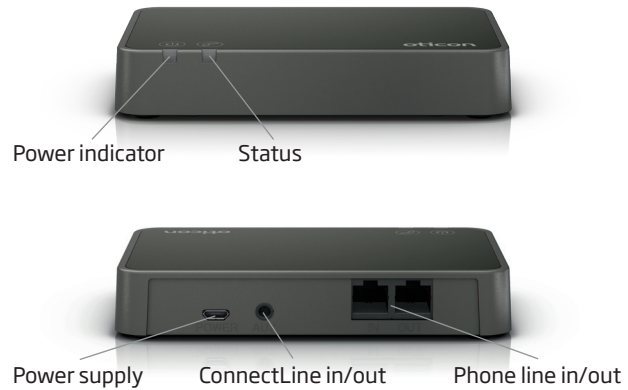

Technical data	
Compatibility	Streamer 1.4 and later
Dimensions	124 x 80 x 21 mm; weight 100 g.
ConnectLine system features	AutoConnect: will automatically connect to the Streamer when in range. AutoPriority: incoming calls will be prioritized over, e.g., TV sound. Proximity pairing: the Phone Adapter will automatically pair to the Streamer if the Streamer is in pairing mode.
Range	Up to 30 meters with Streamer Pro (free line of sight)
Connections	Phone line in/out for interface to (analog) landline home phone (PSTN) (RJ11/RJ11) ConnectLine in/out for direct interface to, e.g., a computer <sup>1</sup> or for interface to the ConnectLine office phone kit <sup>2</sup> (3.5 mm mini jack)
Power consumption	Max. 190 mA (using the included power supply)
Headset in	Max 6dBmV(RMS), 10 K $\Omega$
Headset out	-15 -- -10dBmV(RMS), 10 -100 $\Omega$
Audio Bandwidth	License free 2.4 GHz (ISM band)
AutoConnect (range)	Yes (30 meters)
Audio Quality	Mono / 8 kHz sample frq
Frequency	80 Hz - 3.5 kHz (wireless)
Operating temperature	10 - 40 °C
Storage temperature	-20 - 70 °C
Emissions	Max. 39 dB $\mu$ V/m @ 3 meters. In compliance with all relevant standards

### References:

1. A computer interface cable is included with Streamer Pro and is also available as an accessory.
2. The ConnectLine office phone kit is available as an accessory. It includes a headset interface box and a cable for connection to Phone Adapter 2.0.